

TEST MENU

- Available
- In Development

Tumor

- CEA
- AFP
- CA125
- CA15-3
- CA19-9
- tPSA
- fPSA
- HE4
- CYFRA21-1
- CA72-4
- PG I
- PG II
- G-17
- CA50
- CA242
- SCC
- NSE
- ProGRP

Thrombus

- D-dimer
- TM
- TAT
- FDP
- PIC
- t-PAIC

Anemia

- VB12
- Folate
- Ferritin
- RBC Folate

Fertility

- AMH
- Progesterone
- Testosterone
- Prolactin
- LH
- E2
- FSH
- β -HCG
- Inhibin B
- fTesto
- DHEA-S

Cardiac

- cTnI
- hs cTnI
- CK-MB
- MYO
- BNP
- NT-proBNP
- H-FABP
- ST2
- hs cTnT
- hs-CRP

Inflammation

- PCT
- IL-6
- CRP
- SAA

Infectious Disease

- Dengue NS1 Ag
- Dengue IgG
- Dengue IgM
- HIV Ag/Ab
- HBsAg
- Anti-HCV
- Anti-TP
- H.Pylori (HP)
- HBsAb
- HBeAg
- HBeAb
- HBcAb

Bone Metabolism

- 25-OH Vitamin D Total
- PTH

Immunoglobulins

- IgE
- IgM
- IgG
- IgA

Autoimmune

- ANA
- ASO
- RF
- Anti-CCP
- DsDNA IgG
- Anti-SSA/Ro Antibodies
- Anti-SSB/La Antibodies
- Anti-APS Ab

TORCH

- Toxo IgG
- Toxo IgM
- Rubella IgG
- Rubella IgM
- CMV IgG
- CMV IgM
- HSV-1 IgG
- HSV-1 IgM
- HSV-2 IgG
- HSV-2 IgM

Thyroid

- TSH
- TT3
- TT4
- FT3
- FT4
- Tg
- Anti-Tg
- Anti-TPO
- TRAb

DIABETES

- Insulin
- C-Peptide

Others

- S100 β
- Cortisol
- ACTH

Wondfo

TISENC
Science, Technologies and Future Beyond



Accre 120

AUTOMATIC CHEMILUMINESCENCE IMMUNOASSAY (CLIA) SYSTEM




Wondfo

Guangzhou Wondfo Biotech Co., Ltd.
No.8 Lizhishan Road, Science City,
Huangpu District, 510663, Guangzhou, P.R.China
Tel: (+86) 400-830-8768
Website: en.wondfo.com
E-mail: sales@wondfo.com.cn



© 2025 Guangzhou Wondfo Biotech Co., Ltd.
All rights reserved. Rev.A01 Rel:2025/07/29

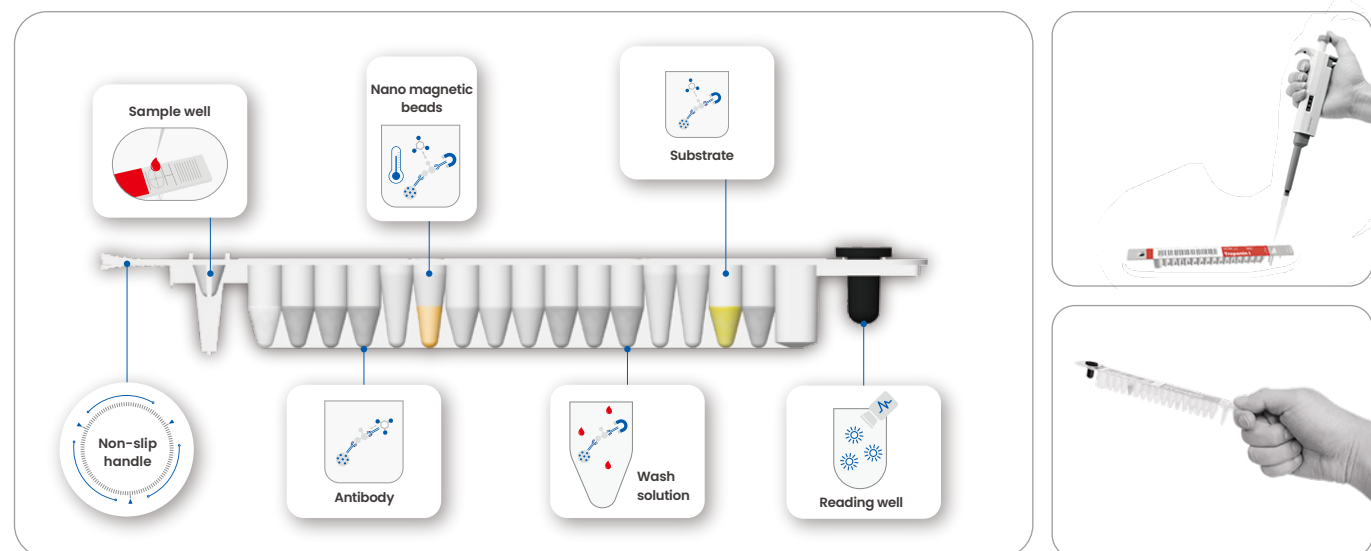

EXCELLENT AS CLIA

- ALP-APS-5 chemiluminescence
- Excellent reagent performance
- Nano magnetic particles for separation
- Broad menu with featured items
- Reliable calibration and quality control system

EASY AS POCT

- Mono-test cartridge, ready-to-use
- Avoid cross-contamination
- Short turnaround time (TAT)
- All-in-one design, no consumables and extra cost
- Built-in scanner reads reagent information automatically

ASSAY PROCEDURE



1 Incubation

Precise temperature
(37°C or 98.6°F) control




2 Separation

Nano-size magnetic particles
for complete separation



3 Reaction

ALP (alkaline phosphatase)-
catalyzed CLIA reaction

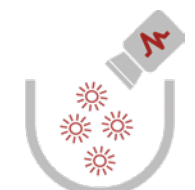

4 Reading

Ultra-sensitive PMT for
signal capture

ACCURATE

CONVENIENT

CONSISTENT

RAPID

ECONOMICAL



Fully Automatic User Experience

- Primary tube sampling



Flexibility & Efficiency

- Simultaneous processing of different parameters
- STAT mode supported
- Higher throughput



Better User Experience

- Automatic barcode identification
- 600 days MTBF

System

Bench-top CLIA analyzer

Methodology

Chemiluminescence Immunoassay (CLIA)

Dimension

570mm*685mm*740mm

Weight

100 kg

Throughput

Up to 120 tests/hour

Result

10 mins for the first result

Mode

Batch/STAT mode

Operation

15.6-inch touch screen, 24 hours standby

Channel

4 modules with 20 channels totally

Sample Type

Whole blood/Serum/Plasma

Maintenance

600 days MTBF^①,
no consumable needed for maintenance

Accessories

Built-in thermal printer and barcode scanner

Data Connection

LIS/HIS connection RS232,
USB, Internet connection

① MTBF: Mean Time Between Failure